




"Buttockmate" type 2

- Intramuscular Injection Simulator

● Weight/About 1.2kg ● Case size/W46×D16×H26cm ● Accessories/Belt, Drain bag, Storage case

Reproducing Clarke's point and Hochstetter's point.
Allowing simulation in the lateral recumbent position.

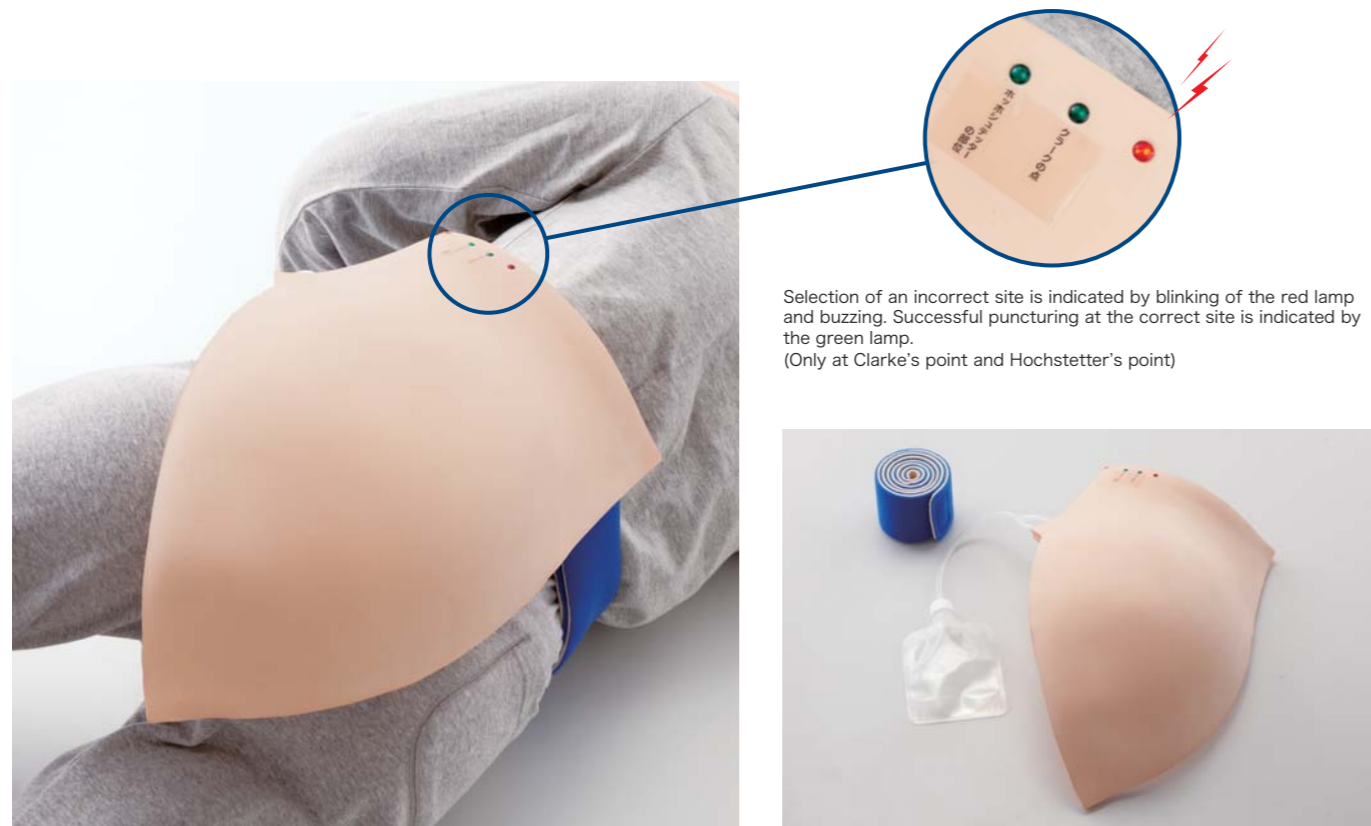

Selection of an incorrect site is indicated by blinking of the red lamp and buzzing. Successful puncturing at the correct site is indicated by the green lamp.
(Only at Clarke's point and Hochstetter's point)

The epidermis is made of a special material that leaves no needle mark.

Being compact, the device needs a small space for storage.

Practice

Hochstetter's point



Skeletal reference points such as the anterior superior and anterior posterior iliac spines and greater trochanter can be felt for efficient training in finding the injection site.

Clarke's point



Fluid can be injected.
* Do not use a drug solution. Use water (tap water only) for practice.



"Buttockmate" - Intramuscular Injection Simulator

● Weight/About 2kg ● Case size/W50×D35×H15cm ● Accessories/Belt, Drain bag, Transparent package, Storage case

Practice gluteus injection skill by palpation. (*Shibusanbu** method)

*Upper Outer quadrant of Buttocks



Feature

- Students can train in pair playing the role of nurse and patient alternately, and this wearing type enables to compare with real body.
- As ilium crista, greater trochanter and iliac spine are built in, confirmation and practice for gluteus injection skill by palpation available.
- Touch of skin and frame that are very close to living body contribute greatly to learn gluteus injection skill.
- The flashing red light and the buzzer warn in case of miss injection to incorrect point or too deep injection though done to correct point. The green light turns on in case of correct point.
- The surface is so similar to the real skin for its specific process that no needle mark remains though repeatedly used.



The drawing of injection site is written on a transparent package.

Practice



By palpation, injection training to correct site is possible.



Water (only tap water) can be infused.



By using the package, practice can be made alone.