


I.V. Trainer with Indicator - Intravenous Injection Trainer

● Weight/About 4kg ● Case size/W51×D42×H20cm ● Accessories/Belt, Exchangeable Blood Vessel, Storage case

Built-in lamps, chime sound and a buzzer indicate correct or incorrect injection into the three main veins: Basilic vein, Cephalic vein and Median cubital vein.

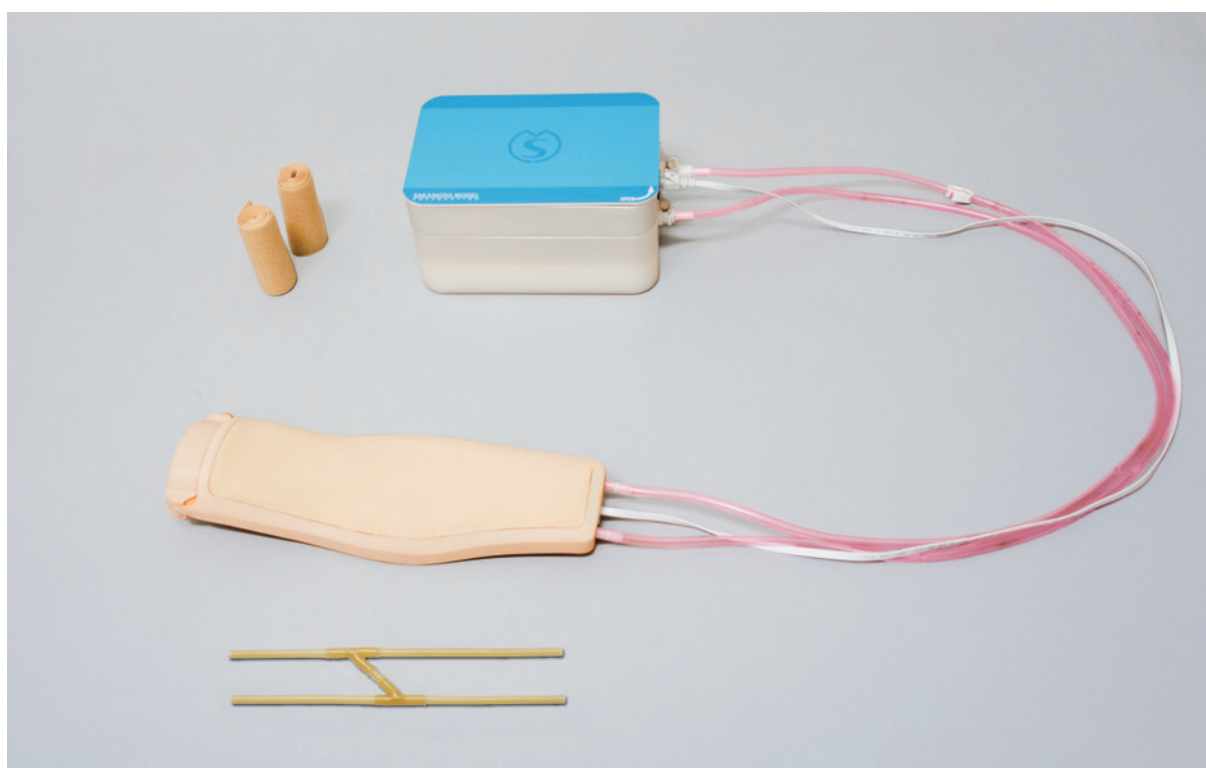

Feature

- 1 Built-in lamps, chime sound and a buzzer indicate correct or incorrect! (Can be used for a test)

OK (Correct injection)

Green lamp, Chime sound




NG (Incorrect injection)

Red lamp, Buzzer



- 2 Blood circulates by an electric pump. Can check back flow.

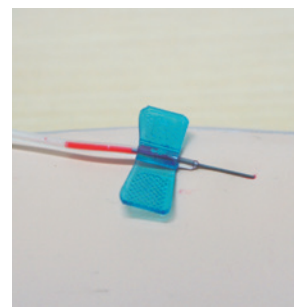

Intravenous injection (puncture)



Blood Sampling



Insertion of indwelling needle/ Fixing with adhesive tape



Exchangeable parts

Easy to change vessels

Exchangeable Blood Vessel



M182

Arm for Simulator

●Size/W60×D25×H23cm

※M182 can be used for M118-1, M148-3, M178, M181, M183, AM-1, AM-2

