

# Sakamoto Acupuncture Simulator for Sciatic Nerve

● Weight/About 3.2kg ● Case size/W60×D47×H41cm

## Practical training using a lateral position model simulating an actual patient!


This model permits the practical training at a lateral position simulating actual treatment.

We aim to learn techniques that allow for accurate insertion using a system that matches the patient's posture.

Since the skin is made of special material, it is possible to learn in a more realistic environment.

The iliac crest (iliac spine), great trochanter, and ischial tuberosity are incorporated in the model. The trainees can practice palpating, and the green light turns on and a chime sounds when the needle is inserted to the correct point facilitating easy self-practice.

Even beginners can practice easily.

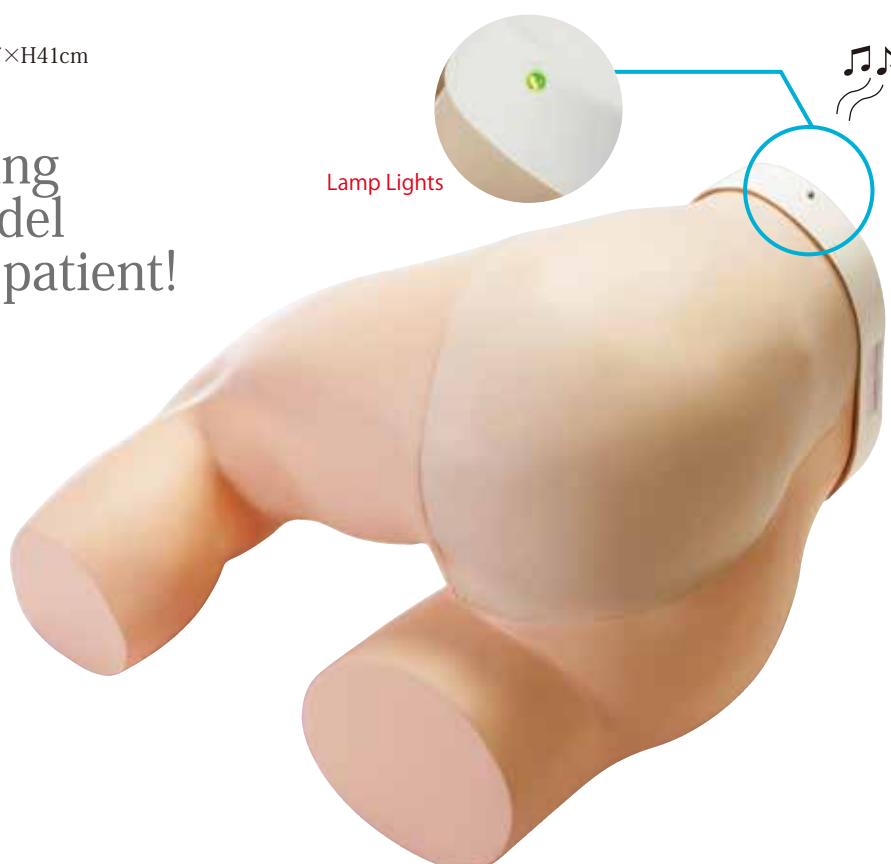

### Practice

#### 1 【Method 1】 Posterior superior iliac spine-greater trochanteric line



Insertion hole: Determine the midpoint between the posterior superior iliac spine and the greater trochanter, and find the insertion point 3 cm below this midpoint at a right angle.

#### 2 【Method 2】 Ichial tuberosity-greater trochanter line



Insertion hole: Connect the greater trochanter and ischial tuberosity, and make a hole at a point 1/3 from the ischial tuberosity side.

#### 3 Acupuncture insertion

You can practice insertion while lying on your side.

