

M156

Decubitus Treatment Simulator

- Quantity/4 types/set ● Weight (1 piece)/About 1.5kg ● Size (1 piece)/W39×D22×H10cm
- Directors/Takaoka station south clinic president Kunio Tsukada, Kanazawa University premedical course Hiromi Sanada

Decubitus series describing particularity of each stage from I~IV with realism

The coloring expresses the characteristics of each stage in a set of 4 types from stages I to IV. Since the coloring is done from the inside, it will not fade even if you use chemicals or wash it. In addition, the pressure wound area does not change color even after training. The structure of epidermal damage and the soft sensation are realistically reproduced for each stage, allowing you to carefully practice care according to each process. It can be attached to the human body or to the model doll with a mounting belt.



Stage I

The epiderm has not yet been torn but a rash can be observed.



Stage II

Abrasion, bleeding, and effusion can be observed.



Stage III

Healing process of stage IV where the fascia and the deeper muscle layer has been affected.



Stage IV

"There exists a hole widely separating the fascia and periosteum from the hypodermal organ."




M195

Suturing Trainer


- Weight: 165g ● Case size: W31×D17×H10cm

The skin can be cut as desired to carry out repeated suturing training

Two types of bases, a flat base and a base which is curved to resemble the human body, can be used for suture training. The skin is replaceable and can be used as many times as you want by using a replacement skin (sold separately).

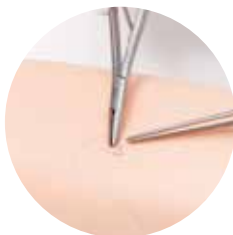

A variety of different suturing training can be carried out.

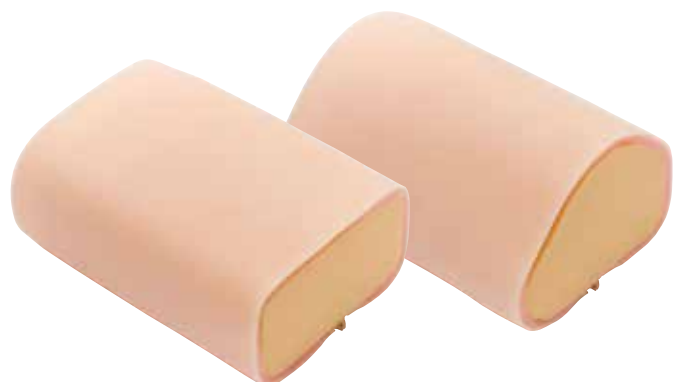

Exchangeable parts

Skin 3 sheet/set